


Original article

Optimal parameters for the enhancement of human osteoblast-like cell proliferation *in vitro* via shear stress induced by high-frequency mechanical vibration

Nahum Rosenberg ^{a,*} , Orit Rosenberg ^b, Jacob Halevi Politch ^c, Haim Abramovich ^c

^a Ruth and Bruce Rappaport Faculty Of Medicine, Technion-Israel Institute of Technology, Haifa, Israel

^b Sheltagen Medical Ltd. Haifa, Israel

^c Department Faculty of Aerospace Engineering, Technion-Israel Institute of Technology, Haifa, Israel

ARTICLE INFO
Article history:

Received 23 March 2021

Received in revised form

21 April 2021

Accepted 05 May 2021

Keywords:

Osteoblast

Mechanical stimulation

Vibration

Mechanotransduction

ABSTRACT

Introduction: Biomechanical stimulation of cultured human osteoblast-like cells, which is based on controlled mechanical vibration, has been previously indicated, but the exact mechanical parameters that are effective for cells' proliferation enhancement are still elusive due to the lack of direct data recordings from the stimulated cells in culture. Therefore, we developed a low friction tunable system that enables recording of a narrow range of mechanical parameters, above the infrasonic spectrum, that applied uniformly to human osteoblast-like cells in monolayer culture, aiming to identify a range of mechanical parameters that are effective to enhance osteoblast proliferation *in vitro*.

Methods: Human osteoblast-like cells in explant monolayer culture samples were exposed to mechanical vibration in the 10-70Hz range of frequencies for two minutes, in four 24 hours intervals. Cell numbers in culture, cellular alkaline phosphatase activity (a marker of cell maturation), and lactate dehydrogenase activity in culture media (representing cell death) were measured after the mechanical stimulation protocol application and compared statistically to the control cell cultures kept in static conditions. The cell proliferation was deduced from cell number in culture and cell death measurements.

Results: We found that 50-70 Hz of vibration frequency protocol (10-30 μ m of maximal displacement amplitude, 0.03g of peak-to-peak acceleration) is optimal for enhancing cells' proliferation ($p < 0.05$), with a parallel decrease of their maturation ($p < 0.01$).

Discussion: We detected the optimal mechanical parameters of excitation protocol for induction of osteoblast proliferation *in vitro* by a mechanical platform, which can be used as a standardized method in the research of mechanotransduction in human osteoblast.

© 2021 The Authors. Published by Iberoamerican Journal of Medicine. This is an open access article under the CC BY license (<http://creativecommons.org/licenses/by/4.0/>).

* Corresponding author.

E-mail address: nahumrosenberg@hotmail.com

ISSN: 2695-5075 / © 2021 The Authors. Published by Iberoamerican Journal of Medicine. This is an open access article under the CC BY license

(<http://creativecommons.org/licenses/by/4.0/>).

<http://doi.org/10.5281/zenodo.4746404>

Parámetros óptimos para la mejora de la proliferación de células similares a osteoblastos humanos *in vitro* a través del esfuerzo cortante inducido por vibraciones mecánicas de alta frecuencia

INFO. ARTÍCULO

Historia del artículo:

Recibido 23 Marzo 2021

Recibido en forma

revisada 21 Abril 2021

Aceptado 05 Mayo 2021

Palabras clave:

Osteoblasto

Estimulación mecánica

Vibración

Transducción mecánica

RESUMEN

Introducción: La estimulación biomecánica de células similares a osteoblastos humanos cultivadas, que se basa en vibraciones mecánicas controladas se ha demostrado anteriormente, pero los parámetros mecánicos exactos que son efectivos para la mejora de la proliferación de células aún son difíciles de alcanzar debido a la falta de registros de datos directos de las células estimuladas en cultivo. Por lo tanto, desarrollamos un sistema sintonizable de baja fricción que permite el registro de un rango estrecho de parámetros mecánicos, por encima del espectro infrasónico, que se aplica de manera uniforme a células similares a osteoblastos humanos en cultivo monocapa, con el objetivo de identificar un rango de parámetros mecánicos que son efectivos para mejorar la proliferación de osteoblastos *in vitro*.

Métodos: Se expusieron células similares a osteoblastos humanos en muestras de cultivo de monocapa de explante a vibración mecánica en el rango de frecuencias de 10-70 Hz durante dos minutos, en cuatro intervalos de 24 horas. El número de células en cultivo, la actividad de la fosfatasa alcalina celular (un marcador de maduración celular) y la actividad de la lactato deshidrogenasa en los medios de cultivo (que representa la muerte celular) se midieron después de la aplicación del protocolo de estimulación mecánica y se compararon estadísticamente con los cultivos de células de control mantenidos en condiciones estáticas. La proliferación celular se dedujo del número de células en cultivo y mediciones de muerte celular.

Resultados: Encontramos que 50-70 Hz de protocolo de frecuencia de vibración (10-30 μm de amplitud de desplazamiento máxima, 0,03 g de aceleración de pico a pico) es óptimo para mejorar la proliferación de células ($p < 0,05$), con una disminución paralela de su maduración ($p < 0,01$).

Discusión: Detectamos los parámetros mecánicos óptimos del protocolo de excitación para la inducción de la proliferación de osteoblastos *in vitro* mediante una plataforma mecánica, que puede utilizarse como método estandarizado en la investigación de la transducción mecánica en osteoblastos humanos.

© 2021 Los Autores. Publicado por Iberoamerican Journal of Medicine. Éste es un artículo en acceso abierto bajo licencia CC BY (<http://creativecommons.org/licenses/by/4.0/>).

HOW TO CITE THIS ARTICLE: Rosenberg N, Rosenberg O, Halevi Politch J, Abramovich H. Optimal parameters for the enhancement of human osteoblast-like cell proliferation *in vitro* via shear stress induced by high-frequency mechanical vibration. Iberoam J Med. 2021;3(3):204-211. doi: 10.5281/zenodo.4746404.

1. INTRODUCTION

Mechanical loading increases bone mass by activating osteoblast metabolic activity and proliferation via cellular biomechanical pathways. These pathways are studied *in vitro* in models of cellular mechanical stimulation [1-3]. The widely used experimental methods for this purpose utilize stretching of cells adherent to elastic membranes or implementation of controlled external fluid flow to cell cultures [1-3]. In both, the mechanism of cellular activation is due to shearing forces and subsequential cellular deformation, which causes cytoskeletal activation, mainly involving the microtubular and microfilament components [4, 5]. These experimental methods require sophisticated hardware with essential special handling and complex experimental protocols, but usually unable to produce uniform stretching forces on all the cells in culture, and

commonly restricted to the mechanical loading frequencies up to 5Hz [6].

Previously a more versatile experimental approach has been proposed, which is highly reproducible and with much wider range of mechanical frequencies that can be applied on the cultured osteoblasts in monolayer. This method is based on the application of controlled vibration force in a defined range of mechanical parameters [7].

In this model, the cellular deformation is caused by shearing forces, applied by extra- an intracellular fluid flow, which is induced by the one-dimensional accelerated vibration movement of cells adherent to a plastic surface. Since the cells in this model are in the same environmental condition, they exposed to the same magnitude of force during the vibration application on a culture well plate. Therefore, large numbers of cells can be studied in the same mechanical conditions, which are determined by the dimensions of the well plate and its weight, its surface

biophysical properties, volume of the culture media in the wells and by a profile of the force generated by a vibration actuator. Additionally, this experimental setup can be used for application of vibration in higher frequencies, i.e., above the infrasonic range (20-70 Hz).

There are several previous reports that human osteoblast *in vitro* can be stimulated by the high-frequency vibration frequencies [4,7,8], probably by mimicking the natural influence of resting muscle on an attached bone, as it is seen on the vibromyogram [9] following muscle's spontaneous periodic contraction in this range of frequencies. These reports are also supported by the *in vivo* studies that showed an enhancing effect of vibration in this range of frequencies on the bone mass buildup [10].

It has been shown that different sub-ranges of effective vibration parameters for osteoblast proliferation and metabolic activity stimulation exist and not overlapping [7]. The main difficulty that evolved in that model was to keep a constant pattern of vibration parameters and to record these values directly from the moving supporting well plate due to mechanical instability of the electromagnetic actuator and the difficulty of the fine-tuning of the mechanical parameters in this range of frequencies. The source of the variability of vibration signal was the surface friction of the moving parts, i.e., in that model, it was impossible to effectively eliminate or reduce the additional vibration frequencies in the higher or lower range, which were produced by the vibration actuator and by the well plate movement. Because the optimal range of vibration frequencies for osteoblast activation is narrow, the delivery of defined stable external vibration signals is crucial. Additionally, to claim that a specific range of vibration frequencies has an enhancing effect on cultured osteoblast, a close to "pure" sinusoidal vibration force pattern should be applied. Therefore, to determine the role of mechanical excitation of osteoblast by high-frequency vibration, we developed an experimental model that utilizes a low friction system to deliver vibration force to a cell culture well plate. This

system provides a relatively stable sinusoidal vibration force that can be effectively controlled and measured.

We hypothesize that by using the designed method of high frequency vibration application, we will be able to provide reliable evidence that human osteoblast is stimulated for enhanced proliferation and will determine the range of on the effective mechanical parameters that trigger this effect.

2. METHODS

2.1. CELL CULTURE PROTOCOL

We studied cultured human osteoblast-like cells primary explant cultures from samples of cancellous bone harvested during knee arthroplasties from distal femurs. The method was approved by the institutional Ethical Committee (No. 1240A). Every donor provided written informed consent. Samples of cancellous bone (2-3 grams each), were incubated in osteogenic media (DMEM with heat-inactivated fetal calf serum (10%), 2mM L-Glutamine, 100µM Ascorbate-2-Phosphate, 20mM HEPES buffer, 10nM Dexamethasone, 150µ/ml Streptomycin ,50 U/ml Penicillin) at 37°C within the humidified atmospheric environment of 95% air with 5% CO₂ (v:v) for 20 - 30 days. The cells grow adherent to the plastic tissue culture plates. This type of cells has specific osteoblast-like characteristics, i.e. expression of alkaline phosphatase, synthesis of osteopontin, polygonal multipolar morphology, positive Von Kossa staining, synthesis extracellular matrix with predominantly type I collagen and small amounts of collagen type III and V, non-collagenous proteins like sialoprotein (BSP) and osteocalcin. Additionally, these cells demonstrate matrix mineralization *in vitro* and bone formation *in vivo* [11-13]. The cells were passage into 24 well-plates. Each well is seeded with 2×10^3 cells.

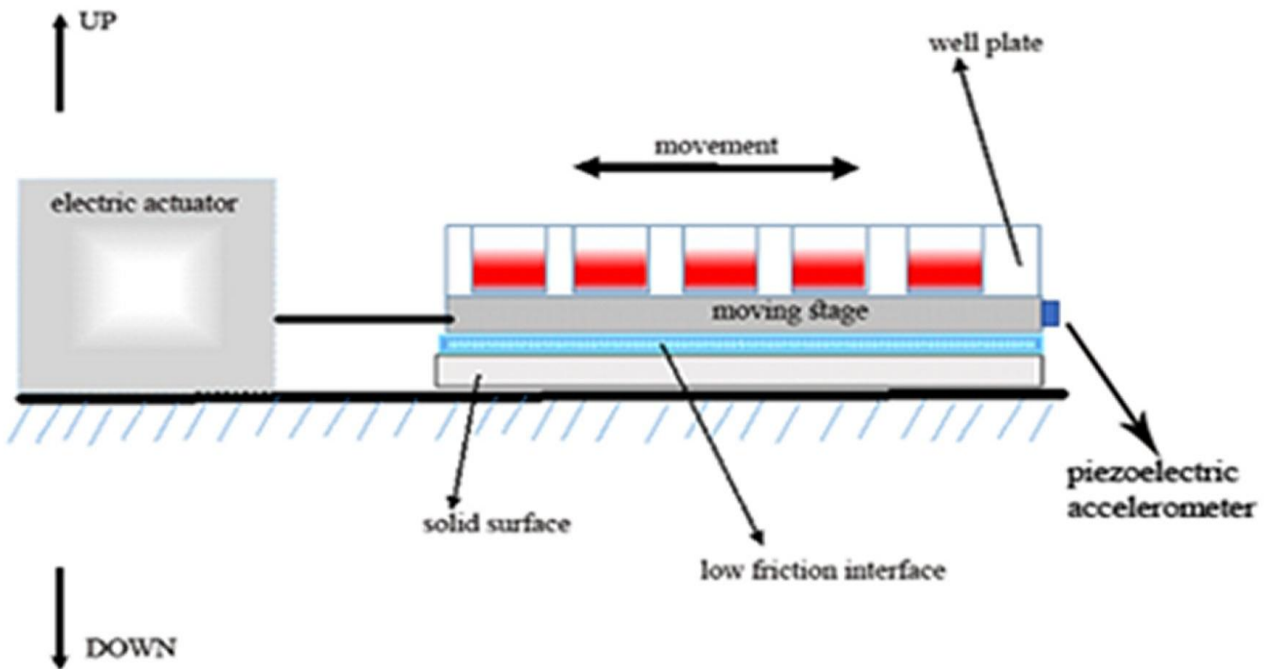


Figure 1: The experimental arrangement (schematic representation). A 24 well-plate containing cells adherent to the plastic surface, covered with culture media, and is tightly attached to the moving stage on a low friction interface. The stage movement is controlled by an electric actuator. The acceleration and displacement of the stage during vibrating movement are recorded by an accelerometer, attached to the stage in a horizontal direction.

2.2. MECHANICAL SYSTEM CONFIGURATION

Plastic 24 well-plate with the cultured osteoblast-like cells (weight 120 grams including the culture media), was mounted on a stage that moves on low friction support (Figure 1). The horizontal movement of the stage with the attached well plate is controlled by an electromagnetic actuator. A piezoelectric accelerometer (accelerometer sensitivity = 100mV/g) is attached to the moving stage with a solid connection to well-plate in the horizontal direction of movement. The relation between the output voltage from a piezo electric accelerometer and the acceleration is linear throughout frequencies of vibration of 10-80 Hz and at the displacement amplitude of up to 30 μ m. This is an open loop, externally controlled setup.

This system produces a vibration movement, which is near sinusoidal, with maximal displacement amplitudes of 25 \pm 5 μ m and a near sinusoidal vibration frequency of 10 and 60 Hz (Figure 2). We defined the dominant and additional harmonic frequencies of the moving well plate from the recordings of the piezoelectric accelerometer (Figure 2). Through a double integration, the displacement of the well plate movement is measured. According to the detected spectra of the vibration applied by this method it was

revealed that when the actuator is tuned to 20 Hz of vibration, 70% of the frequencies are within the range of 10 - 30 Hz (peak to peak acceleration \approx 0.02g), and when the system is tuned to excitation of 60 Hz vibration, 98% of the frequencies are within the range of 50 - 70 Hz (peak to peak acceleration \approx 0.03g).

2.3. MECHANICAL STIMULATION PROTOCOL OF THE CULTURED OSTEOBLAST-LIKE CELLS

Following first passage of cells, replicates of an equal number of cultured cells were submitted to the excitation conditions of 20Hz or 60 Hz. The applied protocol of cultured replicates to mechanical excitation was corresponding to the previous similar studies [4, 7] and utilized two-minute exposure repeated four times at 24 hours intervals. The doubling time of osteoblasts relies on the seeding density. When seeded at the range of 10³ cell/cm² density (as during this study), the doubling time is around 10 hours (14). Therefore, four days of the experiment (96 hours) represent 10 doubling periods. At this time, the cell culture is near 2D confluence, and also the proliferation rate ceases, and the cells start with higher metabolic activity in matrix elaboration. Therefore, within the period of 4 days, the cells reach the maximal proliferation rate. Accordingly, we evaluated the results on

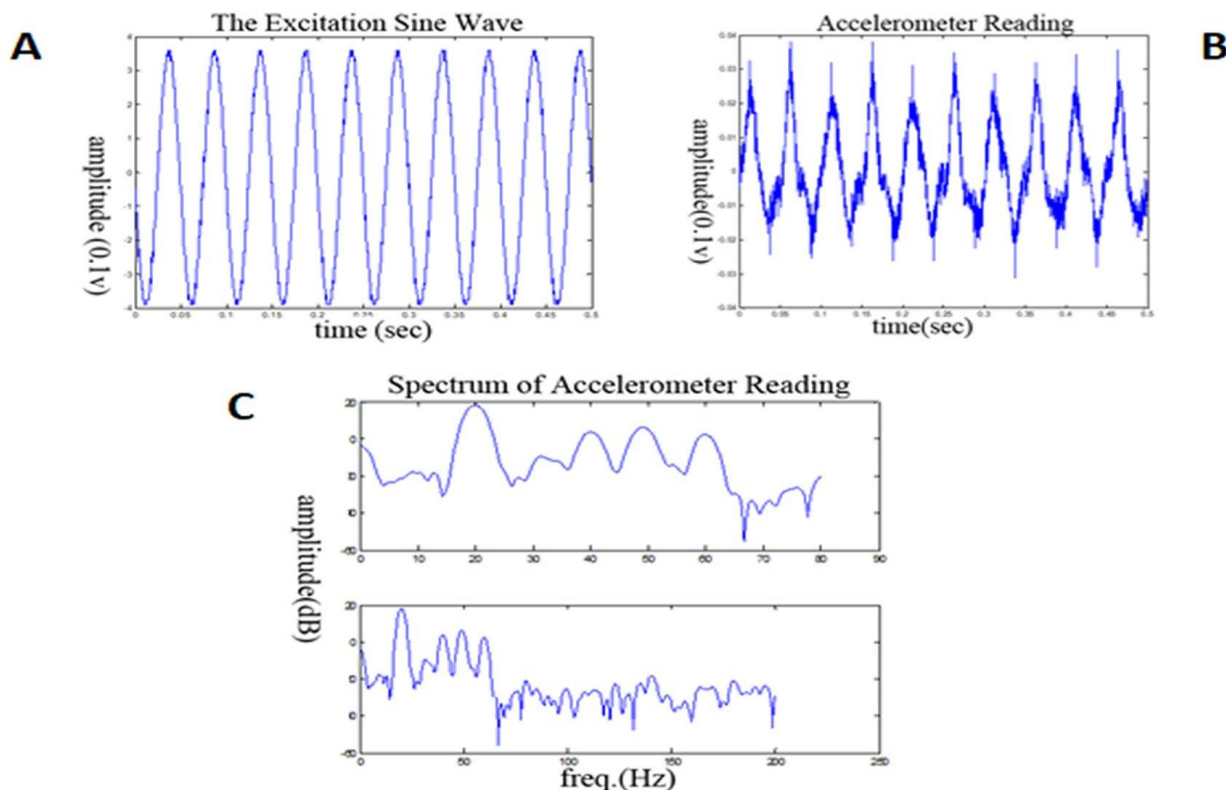


Figure 2: Illustrations of the 20 Hz excitation profile delivered by the actuator (A) and of the accelerometer recording from the moving stage (B). The accelerometer readings reflect the actual mechanical force acting on cells, which are adherent to the moving surface, and appears to be close to similar to the excitation signal from the actuator (amplitude readings represented in $V \times 10^{-1}$). The maximal amplitude of displacement is $25 \pm 5 \mu\text{m}$. The spectra of the vibration frequencies, which were applied to the cells (C) shows that 70% of the frequencies, up to 70 Hz frequency, is in the range of 10-30 Hz.

this point in time of the cell culture growth.

The culture samples kept in static conditions, without exposure to vibration protocols, we used as controls. Eight cell culture samples for every experimental condition were studied.

The number of cells in each culture sample was estimated by direct counting in low-power microscopic field. The average of cell numbers from three different microscopic fields in each culture sample, after exposure to mechanical stimulation and in control samples, was recorded. This method of cell number measurement doesn't require interference with the culture conditions.

2.4. CELL DEATH AND PROLIFERATION

LDH activity within the culture media, a marker of cell death because of LDH leakage via damaged cell membranes, was measured in the samples exposed to mechanical vibration and in the static controls. Briefly, LDH activity in culture media samples was determined by 340 nm wavelength spectrophotometry of the reduced nicotinamide adenine dinucleotide (NAD), which is directly proportional to LDH activity [15, 16].

Accordingly, the cell proliferation comparison among the studied samples in different conditions, i.e., exposed to 60 Hz, 20 Hz of vibration protocols and static controls, was based on as a result of the deduction of cell death change from the total cell number in the culture samples.

2.5. CELLULAR MATURATION

Cellular alkaline phosphatase activity, a marker of osteoblast maturation from progenitors, increases in the Stage 2 osteoblast cell cycle after the Stage 1 cell proliferation ceases [17].

The cellular alkaline phosphatase activity was measured by 410 nm wavelength spectrophotometry [18] after microscopic cell counting, utilizing cell lysate following incubation with P-nitrophenyl phosphate substrate.

2.6. DATA ANALYSIS

The comparison of two parameters was made by the *t*-test and the comparison of three parameters by ANOVA, and subsequently by the Dunnett's post hoc test (following verification of normal distribution of results and with a *p*-

value below 0.05). A P-value of less than 0.05 was considered as a significant difference between the compared values. The calculations of the statistical comparisons were done using SigmaStat software (version 2, SPSS Inc., Chicago, IL, USA).

3. RESULTS

At the 60 Hz vibration the number of cells significantly increased in comparison to control static condition (average 68 ± 4 SEM cells/microscopic field vs. average 55 ± 4 SEM cells/microscopic field respectively, $p < 0.05$, one-way ANOVA, $n=8$, Figures 3 and 4). No significant change in number of cells exposed to 20Hz protocol in comparison to static controls was observed (average 64 ± 3 SEM cells/microscopic field vs. average 55 ± 4 SEM cells/microscopic field respectively, $p > 0.05$, one way ANOVA, $n(\text{controls})=8$, $n(20\text{ Hz})=12$, Figures 3 and 4).

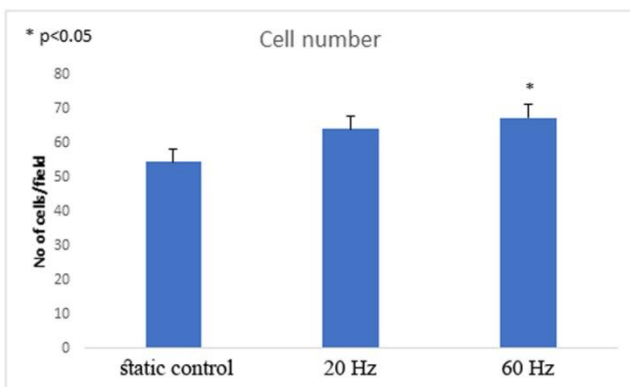


Figure 3: The number of cells rises significantly following 60 Hz excitation (average 68 ± 4 SEM cells/microscopic field vs. average 55 ± 4 SEM cells/microscopic field in controls, $p < 0.05$, $n=12$).

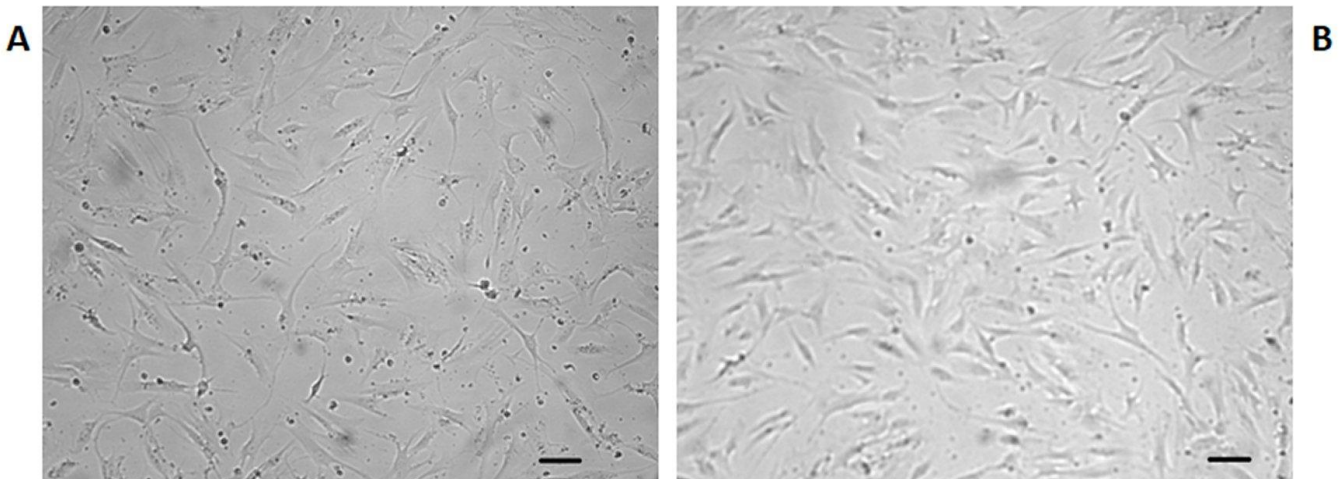


Figure 4: Representative micrographs of osteoblasts in culture on the fourth day of the experiment (Scale 20 μm): A: Cells kept in static conditions. B: Cells after exposure to 60 Hz vibration protocol. More cells in B are apparent.

There was no difference in the media LDH activity in cultures exposed to either 20Hz or 60Hz vibration protocols and the static controls (average 163 ± 9 SEM U/L and average 167 ± 8 SEM U/L vs. average 193 ± 12 SEM U/L respectively, $p > 0.05$, one-way ANOVA, $n=12$).

A significant decrease in cellular alkaline phosphatase activity following 60 Hz vibration protocol was found (average 0.559 ± 0.041 SEM U/L/cell in the microscopic field vs. average 0.794 ± 0.0571 SEM U/L/cell in microscopic field in static controls, $p < 0.01$, t-test, $n=12$, Figure 5).

4. DISCUSSION

In this study, we determined the narrow range of high frequencies spectrum of 50-70 Hz mechanical vibration (when excitation by 60Hz provided by the actuator) that is effective in enhancing human osteoblast proliferation in vitro. This conclusion is following the deduction of cell death rate from cell numbers per sample in different mechanical protocols that were studied, i.e., 60 Hz excitation of cells caused a significant increase in cell number in comparison to controls, without a change in cell death (expressed by LDH activity in culture media) meaning that an enhancement in cell proliferation has occurred. This phenomenon was not apparent following 20Hz excitation. Furthermore, this conclusion is independently supported by the evidence of the lower maturation state of the cells exposed to the 50-70 Hz mechanical vibration protocol as expressed by lower alkaline phosphatase cellular activity. We preferred this double approach of independent assays to increase the confidence in the results. The results show that after the exposure to the around 60 Hz excitation, the cells react by

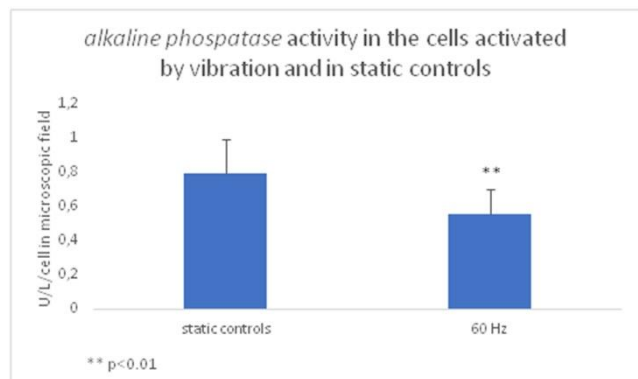


Figure 5: After the 60 Hz excitation, the cells' maturation state decreased -significantly lower alkaline phosphatase activity in the mechanically activated cells was found (average 0.559 ± 0.041 SEM U/L/cell in the microscopic field vs. average 0.794 ± 0.0571 SEM U/L/cell in microscopic field, $p < 0.01$, $n = 12$).

increased proliferation and stay in a lower differentiation state, as expected. This effect doesn't happen in the lower range of the vibration parameters (cellular excitation by 20 Hz that provide 10-30 Hz of vibration spectrum), i.e., by this excitation protocol neither cell number no LDH activity in culture media were not significantly different from the cells kept in control static environment. This fact supports our basic hypothesis that above 10 Hz of mechanical cellular excitation only a specific band of frequencies of vibration elicits a proliferation response. This observation is different from the previous report on a more pronounced proliferation response by 20Hz excitation of human osteoblasts [7], but similar to a previous report on enhanced osteoblast proliferation following 45 Hz vibration protocol [19]. We attribute this controversy to the differences in the experimental setups in these previous reports in terms of acceleration magnitude and direction, time of the vibration protocols applications, and the methods of mechanical measurements. In the present study we determined the narrow range of mechanical frequency of high-frequency horizontal vibration, which causes the highest rate of cell proliferation in human osteoblast that has not been defined precisely till now.

The mechanotransduction in osteoblast involves numerous cellular pathways that propagate mechanical force from the extracellular milieu and translated to cellular synthetic response. The exact interrelation of these pathways is not sufficiently clear and is under extensive research efforts. Overall, this process is initiated by the external mechanical force causing shear stress on the cellular membrane and governed by the mechanical stiffness properties of the supporting extracellular matrix. The deformation of the cellular membrane affects the cell membrane potential by inducing electrical currents of extremely low

electromagnetic intensities (0.02-0.03 mT) via cellular membrane by a flux of calcium and phosphate ions [20-22]. Additionally, the external force is transferred via transmembrane molecular extensions (Integrins, Cadherins, etc.) and cytoskeletal components to a nuclear Lamin A with a subsequential transcription activity [23].

In this study, we did not attempt to determine the specific mechanotransduction biomechanical pathways that are evoked by this mechanical excitation method, but rather provide the basic optimal biomechanical data for the further mechanotransduction research by the presented experimental setup that generates a uniform force transfer to a large number of cells in two-dimensional culture.

The main limitation of this study is the lack of a determination of the possible change of the cell cycle profile following the discovered enhancing proliferative effect of 50-70 Hz mechanical vibration. This important issue can substantially support the presented here experimental results and can be further investigated by the presented experimental setup, which is easily reproducible.

5. CONCLUSIONS

According to the study hypothesis, we describe a tunable and easily reproducible mechanical platform that can be utilized for the research of mechanotransduction in osteoblasts, with the determination of the optimal mechanical frequencies for induction of cell proliferation. Naturally, other cellular functions related to mechanical stimulation can be studied by using this versatile method.

6. REFERENCES

1. Jung S, Bohner L, Spindler K, Hanisch M, Kleinheinz J, Sielker S. Mechanical stimulation increases the proliferation and differentiation potential of human adipose-derived stromal cells. *Int J Stem Cell Res Ther.* 2018;5(2): 1-7. doi: 10.23937/2469-570X/1410056.
2. Schaffer JL, Rizen M, L'Italien GJ, Benbrahim A, Megerman J, Gerstenfeld LC, Gray ML. Device for the application of a dynamic biaxially uniform and isotropic strain to a flexible cell culture membrane. *J Orthop Res.* 1994;12(5):709-19. doi: 10.1002/jor.1100120514.
3. Soejima K, Klein-Nulend J, Semeins CM, Burger EH. Calvarial and long bone cells derived from adult mice respond similarly to pulsating fluid flow with rapid nitric oxide production. *Calcif Tissue Int.* 1999; 64:S113
4. Rosenberg N. The role of the cytoskeleton in mechanotransduction in human osteoblast-like cells. *Hum Exp Toxicol.* 2003;22(5):271-4. doi: 10.1191/0960327103ht3620a.
5. Toma CD, Ashkar S, Gray ML, Schaffer JL, Gerstenfeld LC. Signal transduction of mechanical stimuli is dependent on microfilament integrity: identification of osteopontin as a mechanically induced gene in osteoblasts. *J Bone Miner Res.* 1997;12(10):1626-36. doi: 10.1359/jbmr.1997.12.10.1626.
6. Steward AJ, Cole JH, Ligler FS, Loba EG. Mechanical and Vascular Cues Synergistically Enhance Osteogenesis in Human Mesenchymal Stem Cells. *Tissue Eng Part A.* 2016;22(15-16):997-1005. doi: 10.1089/ten.TEA.2015.0533.
7. Rosenberg N, Levy M, Francis M. Experimental model for stimulation of cultured human osteoblast-like cells by high frequency vibration. *Cytotechnology.* 2002;39(3):125-30. doi: 10.1023/A:1023925230651.
8. Wang L, Hsu H, Xian J. Vibration analysis of a single osteoblast in vitro using the finite element method. *Int Conf Biomed Biol Eng.* 2016;380-5. doi: 10.2991/bbe-16.2016.58.
9. Nigg BM. Acceleration. In: Nigg NM, Herzog W, editors. *Biomechanics of the musculo-skeletal system.* 2nd ed. Chichester: John Wiley & Sons; 1998: 300-1.
10. Rubin C, Li C, Sun Y, Fritton C, McLeod K. Non-invasive stimulation of trabecular bone formation via low magnitude, high frequency strain. *41st Orthop Res Soc.* 1995;20:548.
11. Rosenberg N, Soudry M, Rosenberg O, Blumenfeld I, Blumenfeld Z. The role of activin A in the human osteoblast cell cycle: a preliminary experimental in vitro study. *Exp Clin Endocrinol Diabetes.* 2010;118(10):708-12. doi: 10.1055/s-0030-1249007.
12. Gundle R, Stewart K, Screen J, Beresford JN. Isolation and culture of human bone-derived cells. In: Beresford N, Owen ME, editors. *Marrow stromal cell culture.* Cambridge, UK: Cambridge University Press; 1998:43-66.
13. Yamanouchi K, Satomura K, Gotoh Y, Kitaoka E, Tobiume S, Kume K, et al. Bone formation by transplanted human osteoblasts cultured within collagen sponge with dexamethasone in vitro. *J Bone Miner Res.* 2001;16(5):857-67. doi: 10.1359/jbmr.2001.16.5.857.
14. Megat Abdul Wahab R, Mohamed Rozali NA, Senafi S, Zainol Abidin IZ, Zainal Ariffin Z, Zainal Ariffin SH. Impact of isolation method on doubling time and the quality of chondrocyte and osteoblast differentiated from murine dental pulp stem cells. *PeerJ.* 2017;5:e3180. doi: 10.7717/peerj.3180.
15. Gay RJ, McComb RB, Bowers GN Jr. Optimum reaction conditions for human lactate dehydrogenase isoenzymes as they affect total lactate dehydrogenase activity. *Clin Chem.* 1968;14(8):740-53.
16. Rosenberg N, Hamoud K, Rosenberg O. Quantitative expression of cell death by LDH activity. *IOSR J Pharm Biol Sci.* 2016;11:46-8.
17. Rutkovskiy A, Stensløkken KO, Vaage IJ. Osteoblast Differentiation at a Glance. *Med Sci Monit Basic Res.* 2016;22:95-106. doi: 10.12659/msmbr.901142.
18. Bessey OA, Lowry OH, Brock MJ. A method for the rapid determination of alkaline phosphates with five cubic millimeters of serum. *J Biol Chem.* 1946;164:321-9.
19. Gao H, Zhai M, Wang P, Zhang X, Cai J, Chen X, et al. Low-level mechanical vibration enhances osteoblastogenesis via a canonical Wnt signaling-associated mechanism. *Mol Med Rep.* 2017;16(1):317-24. doi: 10.3892/mmr.2017.6608.
20. Yost MG, Liburdy RP. Time-varying and static magnetic fields act in combination to alter calcium signal transduction in the lymphocyte. *FEBS Lett.* 1992;296(2):117-22. doi: 10.1016/0014-5793(92)80361-j.
21. Liburdy RP. Calcium signaling in lymphocytes and ELF fields: evidence for an electric field metric and a site of interaction involving the calcium ion channel. *FEBS Lett.* 2000;478(3):304. doi: 10.1016/s0014-5793(00)01541-6.
22. Sundelacruz S, Moody AT, Levin M, Kaplan DL. Membrane Potential Depolarization Alters Calcium Flux and Phosphate Signaling During Osteogenic Differentiation of Human Mesenchymal Stem Cells. *Bioelectricity.* 2019;1(1):56-66. doi: 10.1089/bioe.2018.0005.
23. Khilan AA, Al-Maslmani NA, Horn HF. Cell stretchers and the LINC complex in mechanotransduction. *Arch Biochem Biophys.* 2021;702:108829. doi: 10.1016/j.abb.2021.108829.